


DosiPump DP1000

Intelligent Controlled
Peristaltic Pumps

Touch'n'Pump – easy programming and low pulsation

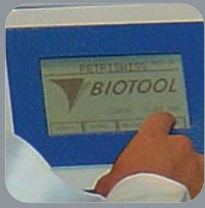

Data Input



Easy - Load



10 Roller Head



Low Pulsation



Low Pulsation

Incorporating a unique 10-roller pump head, the DosiPump eliminates pulsation effects that are incumbent with many conventional peristaltic pumping systems. The result is a precise and reproducible liquid dose for volumes up to 10 litres with a +/- 0.1 ml tolerance.

Versatile and Easy to Use

Capable of single-volume dosing, dilutions, pipetting and multiple dosing. The easy-to-operate DosiPump provides the versatility to undertake any application requiring fast, reproducible and precise dispensing of liquids. Equipped with intuitive and intelligent software. Accessible via a convenient touch screen interface, any operator can be productively using the DosiPump after minimal instruction.

Sturdy Housing

Affordably priced the high performance DosiPump incorporates AntiDrop drip suppression technology preventing dripping between dispensing operations thus reducing the risk of contamination. Built to the highest Swiss quality standards, the DosiPump is simple to maintain and easy to clean and resistant to chemicals.